

MARIOTTI

UNIKO III
Surgery

IMPLANTOLOGY - SURGERY MOTORS

+TORQUE/SPEED PERFORMANCE + MORE PRACTICALITY + MORE RELIABILITY

The improved versions MiniUniko, designed and manufactured in Italy, combine ease of use and practicality with top torque performance, full safety and operating precision with any implant system.

MUN.CL

LED
high luminosity light



MUN.C CLASSIC



MUN.F EASY



HANDPIECES


C20L optic contra-angle for led motors: 20:1 reduction, external and internal irrigation. Top reliability and high torque.



C20 classic contra-angle: 20:1 reduction, external and internal irrigation. Top reliability and high torque.

Constant top performance

- › Newest generation brushless motor and electronic control
- › Accurate adjustment of speed and torque up to 40,000 rpm and 80 Ncm
- › Extended service life of parts that require autoclave sterilization: motor, cable and connector

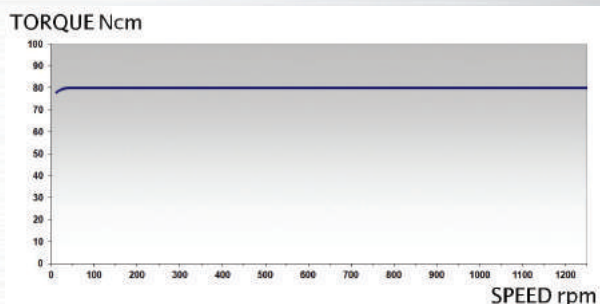

new generation brushless motor with more power




special protection for sterilization

Rev. 02 - 2016

Immediate use and practicality

- › New practical PERISTALTIC PUMP: fast tube replacement and efficiency
- › Easy touch keyboard and wide display
- › Multifunction foot-control with accurate variable speed pedal, reverse, irrigation pump, programs change. (optional for MUN.F version)

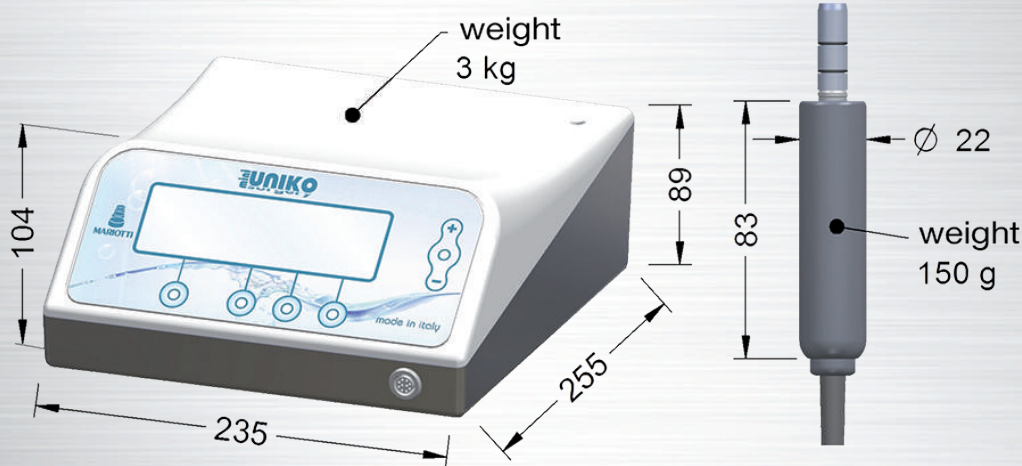

Main Specifications:

- Power supply: 230 V - 115 V ; 50/60Hz
- Absorbed power: 109 VA
- Precise Torque adjustment: max 80Ncm (32:1), 70 Ncm (20:1)
- Precise Speed adjustment: 400-40000 rpm
- Insulation: Class I, type BF
- Peristaltic pump max irrigation: 90 mL/min.
- 10 settable programs
- Irrigation and led functions, automatically/separately activated with the motor start

The High Quality Made in Italy



MARIOTTI & C. SRL
 Implant-surgery motors and piezoelectric - manufacturer since 1973
 Via Seganti 73 - 47121 Forlì - Italy
 Tel +39 0543 474105 - Fax +39 0543 781811
 info@mariotti-italy.com
 www.mariotti-italy.com



T +34 93 419 47 70
 info@adin-iberica.com
 www.adin-iberica.com



Distribuido por:

ADIN IBÉRICA IMPLANTES
 Via Manuel de Falla, Nº 14
 08753 - Fontpineda - Spain